

Home

Welcome to the documentation for the Bonsai workflow developed for the Sound Lateralization Task that is going to be performed by the Circuit Dynamics and Computation group at the Champalimaud Foundation.

The experimental setup makes use of the capabilities of the Harp devices (which implement the [Harp](#) protocol) and the [Bonsai](#) visual programming framework, which work really well together.

The documentation for the task is divided in 3 sections:

- [Installation and Configuration](#) - This section is a step-by-step guide for setting up everything that is needed to run the task in a new setup. It includes instructions for installing the task's software, the firmware for each Harp device used in the task, calibrating the hardware and configuring the task.
- [Task](#) - This is the section where the task is described. It contains a high-level explanation of the task, as well as a more low-level one.
- [Extensions](#) - Since there was a need to create custom Bonsai nodes written in C# to implement functionality which was difficult to implement or wasn't available natively in Bonsai (namely the reading of configuration files), a section documenting these nodes had to be created. There is also a description of every parameter from each configuration file.

If you find any bug in the project or any missing/incorrect/out-of-date documentation, feel free to create an issue on [GitHub](#) or contribute with a pull-request.

If you want to build the documentation locally click [here](#).

Bill of Materials

This bill of materials (BOM) contains **EVERYTHING** needed to build a setup. From the less memorable screw to the super expensive syringes used for reward delivery.

NOTE

This BOM is destined to the people from the Champalimaud Foundation (CF). Nonetheless, the goal is that people outside CF are also able to order everything and build the setup.

WARNING

This BOM specifies the part numbers used by the lab for these setups, but not every component needs to be of a specific part number.

The only things that need to be of a specific model are the Harp Devices, the Behavior Poke Port Breakout v1.1, the Poke Small v1.1 and the 12V Power Supplies.

It's possible to use a different camera model, but please confirm that the camera can be used in Bonsai (although for some cameras the Bonsai workflow will probably need to be adapted).

Harp

This section contains the setup components that relate to the Harp devices somehow (except for the Camera that has a section of its own).

Devices

The Harp devices are assembled by the Hardware and Software Platform, so place an Agendo request to order them.

WARNING

For people outside of the Champalimaud Foundation, it's possible to order the Harp devices in the [Open Ephys Production Site](#)[↗].

It's also possible to assemble the devices in-house with the right equipment since all of them are under an open source license. Some devices (like the Harp SyringePump) must be built in-house since they are not sold externally.

Item	Description	Amount	Observations
Harp Behavior	General-purpose Harp board	1	-
Harp SoundCard	Delivers the auditory stimulus	1	-
Harp Audio Amplifiers	Amplifies the auditory stimulus	2	1 Harp Audio Amplifier per speaker
Harp ClockSynchronizer	Synchronizes the timestamps from every Harp device	1	-
Harp SyringePump	Device used for reward delivery	2	-
Harp CurrentDriver	Drives/controls the LED/laser used in optogenetics	1	Optional



Peripherals

For this part of the BOM, a peripheral is considered to be anything that either interacts with the Harp devices or is needed make the different devices work and/or interact with each other.

Item	Description	Amount	Part Number	Part of the Harp Device kit	Agendo	Observations
Behavior Poke Port Breakout v1.1	Makes the ethernet ports pins from the Harp Behavior available	3	-	✗	✓	Connects to the Harp Behavior
RJ-to-RJ cables	-	3	IM1037	✗	✗	Connects the Harp Behavior to the Behavior Poke Port Breakout v1.1
Poke Small v1.1	Board with infrared beam to detect animal pokes	3	-	✗	✓	Connects to the Behavior Poke Port Breakout v1.1

Item	Description	Amount	Part Number	Part of the Harp Device kit	Agendo	Observations
3.5 mm-stereo-audio-jack-to-bare-wires	-	3	BC-A3ML006F	✗	✗	Connects the Poke Small v1.1 to the Behavior Poke Port Breakout v1.1
5 mm white/blue LED	Placed on the box lid to give cues to the animal	1	-	✗	✗	Connects to the LED0 pins of the Harp Behavior
3 mm green LED	Placed in the central poke to give cues to the animal	1	-	✗	✗	Connects to the LED1 pins of the Harp Behavior
RCA-to-RCA cables	-	2	-	✓	✗	Connects each Harp Audio Amplifier to the Harp SoundCard
Speakers	Deliver the auditory stimulus	2	TODO	✗	✗	1 speaker per Harp Audio Amplifier
Banana Plug	4 mm Triple Contact Plug (Black or Red) - Pack of 10 units	4	557-0100	✗	✗	Connects the speakers to the Harp Audio Amplifiers
10 ml Glass Syringe	Hamilton 1000 Series Gastight Syringes: Luer Lock Syringes,	2	Hamilton 81620	✗	✗	Used for the Harp SyringePump

Item	Description	Amount	Part Number	Part of the Harp Device kit	Agendo	Observations
	TLL Termination					
4-way Stopcock 	Pack of 10 units	2	Masterflex 30600-04	✗	✗	Attaches to the end of the syringe (1 per Harp SyringePump)
Nylon Male Luer Fitting 	Pack of 25 units	2	Masterflex MFLX45505-31	✗	✗	Attaches to one end of the stopcock (1 per Harp SyringePump)
Nylon Female Luer Fitting 	Pack of 25 units	2	Masterflex 45502-00	✗	✗	Attaches to one end of the stopcock (1 per Harp SyringePump)
Spouts 	Used in reward delivery	2	B14200 100 BULK	✗	✗	Used in reward delivery. Glued to the physical lateral pokes
Flexible Tubing 	Masterflex Tygon E-3603 Non-DEHP Tubing - 15 meters per unit	2	Masterflex 06407-71	✗	✗	Connects the spouts to one of the luer fittings (either male or female)
Mini USB cable 	-	5/6*	AK-300130-018-S	✓	✗	Connects the Harp boards to the computer
Micro USB cable	-	1	-	✓	✗	Used to upload sounds to the Harp SoundCard
3.5 mm-stereo-	-	4/5*	TODO	✗	✗	Connects the Harp ClockSynchronizer to

Item	Description	Amount	Part Number	Part of the Harp Device kit	Agendo	Observations
audio-jack-to-jack cable						every other Harp device
12V Power Supply  **	AC/DC Wall Mount Adapter 12V 12W	9/10*	VER12US120-JA	✓ ***	✗	1 is used to power the LED strip that illuminates the behavior box, the remaining are used to power the Harp devices
BNC-to-bare-wires connector  **		1	810-4605	✗	✗	Connects the LED strip to the power supply

* if the Harp CurrentDriver is being used

** the number of power supplies can be decreased as explained in the subsection [below](#)

*** the power supply used to power the LED strip that illuminates the behavior box must be ordered separately

Reducing the number of power supplies

From the table above, it's pretty noticeable that a lot of power supplies are required, which takes up a lot of space and power outlets. It's possible to use the same power supply for different devices according to their characteristics. The devices can be grouped in the following way:

- Harp Behavior, Harp SoundCard and Harp CurrentDriver*
- V+ of the Harp Audio Amplifiers
- V- of the Harp Audio Amplifiers
- Both Harp SyringePumps

From the list above, it's possible to switch the last line from the previous table with the lines from the table below.

Item	Description	Amount	Part Number	Comes with Harp Device	Agendo	Observations
12V Power Supply	AC/DC Wall Mount Adapter 12V 12W	5	VER12US120-JA	✓***	✗	1 is used to power the LED strip that illuminates the behavior box, the remaining are used to power the Harp devices
2-to-1 cables for the 12V power supplies	-	4/3*		✗	✓***	-
4-to-1 cables for the 12V power supplies	-	0/1*		✗	✓***	-

* if the Harp CurrentDriver is being used

Mechanical Components

This section contains the mechanical components of the setup that were developed and assembled in-house.






⚠ WARNING

The files for the mechanical components are not currently available online and can't also be ordered, so people outside of the Champalimaud Foundation will have to develop and assemble their own.

Item	Description	Amount	Part Number	Observations
Behavior box	-	1	-	Made of acrylic
Physical pokes	-	3	-	Preferably made of metal, but can also be 3D printed
Speaker holder	-	2	-	Preferably made of metal, but can also be 3D printed
Speaker holder pole	Hollow aluminium tube	2	-	-
Box LED holder	-	1	-	3D printed

Camera

This section contains the hardware needed to setup the camera and fix it to the lid of the behavior box. Click [here](#) to go to the camera configuration instructions.


Item	Description	Amount	Part Number	Observations
FLIR Camera 		1	BFS-U3-16S2M-CS	-
Camera Lens 		1	A4Z2812CS-MPIR	-
Camera USB cable 	USB-A to Micro-B Straight with Thumbscrews, 3 Meters	1	CEI USB3-1-1-2-3M	This cable MUST be connected to a USB 3.0 port for performance
Camera GPIO cable 	6 Pin Female Straight Plug (Hirose HR10A-7P-6S) to Flying Leads, 3 Meters	1	CEI MVA-50-3-X-3	Connects to the Harp Behavior to trigger/monitor the camera frames
Tripod Adapter 	BFS 30 mm BFLY CM3 Tripod Adapter	1	ACC-01-0003	Used to fix the camera to the Thorlabs poles in the behavior box
Thorlabs 75mm post	Ø12.7 mm Optical Post, SS, M4 Setscrew, M6 Tap, L = 75	1	TR75/M	Connects to the camera

Item	Description	Amount	Part Number	Observations
✚	mm			
Thorlabs 150mm post ✚	Ø12.7 mm Optical Post, SS, M4 Setscrew, M6 Tap, L = 150 mm	1	TR150/M	Fixes the camera structure to the box lid
Thorlabs post clamp ✚	Rotating Clamp for Ø1/2" Posts, 360° Continuously Adjustable, 5 mm Hex	1	SWC/M	Fixes both Thorlabs posts to each other
M6 Screw	15/16 mm	1		Fixes the longer post to the behavior box lid
M6 Setscrew ✚	M6 x 1.0 Stainless Steel Setscrew, 16 mm Long, 25 Pack	1	SS6MS16	Connects the Thorlabs post to the Camera

Others


This section contains the remaining components needed for the setup. None of these components needs to be the exact model present on the list. It just corresponds to the models that have been used by the lab, but feel free to use different ones, as long as they work.

Item	Description	Amount	Part Number	Observations
Computer ✚	Mini PC Blackview MP100 Mini PC AMD Ryzen 7 5825U/16GB/512GB SSD	1	MP100(16+512)-BLACK	-
Screen	-	1	-	-
HDMI cable	-	1	-	-
Keyboard	-	1	-	-
Mouse	-	1	-	-
USB Hub ✚	HUB USB 3.0 TP-Link UH700 7 Ports	1	P004616	-
Power extension with 6 electrical	Must be appropriate to 90° plugs	1		To plug all of the power supplies in a setup

Item	Description	Amount	Part Number	Observations
outlets				
KVM switch 		-		Optional: useful in case one wants to use the same screen + mouse + keyboard kit in different computers
12V LED Strip		1		-
Optogenetics Light Source	-	-	-	Optional

Introduction

The sound lateralization task implemented in the current project, which is based on [1], was designed as a state machine, where the progression through the different states is driven by certain events. The figure below is a representation of the state machine that describes this task.

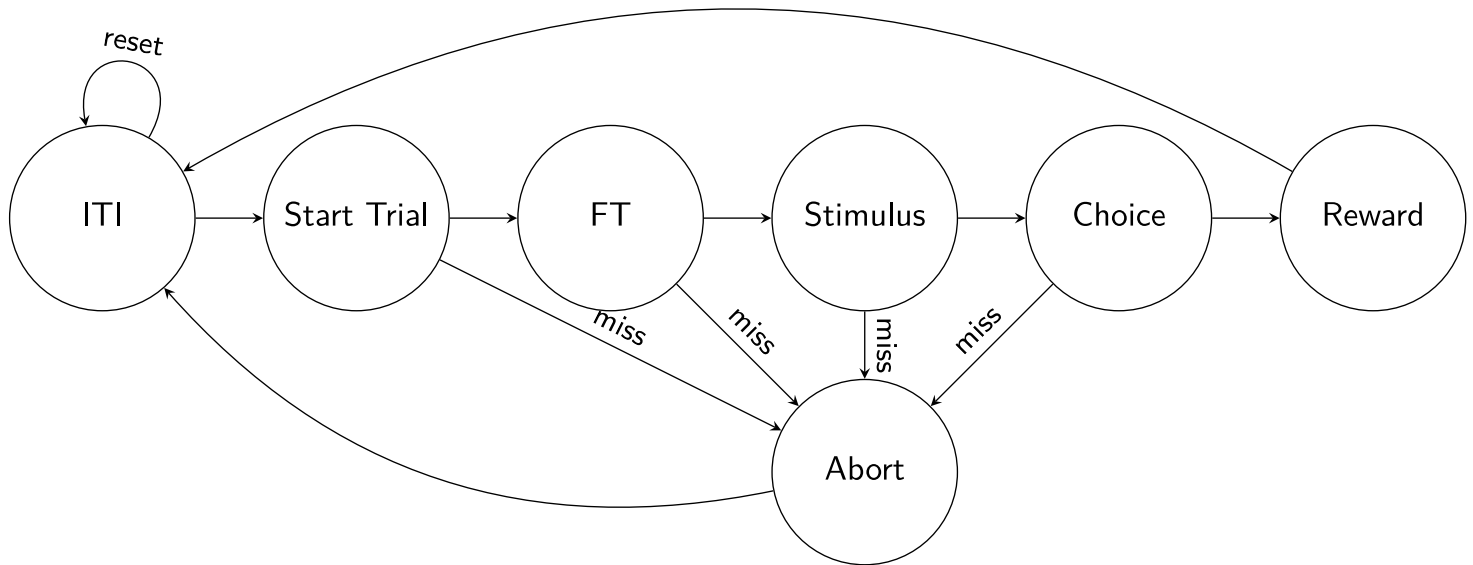

From the figure, notice that from most states there are two possible states that these states can progress to. This happens because there are certain conditions that have to be met in order for the state machine to progress to the next "desired" state.

The remaining pages of this subsection are dedicated to the description each individual state (which include an explanation of what happens and what are the progression conditions).

Namespace Animal

Classes

[Animal](#)

[AutobiasCorrection](#)

[DeserializeFromJson](#)

Deserializes a sequence of JSON strings into data model objects.

[DeserializeFromYaml](#)

Deserializes a sequence of YAML strings into data model objects.

[FixationTime](#)

[OptoLED](#)

[Optogenetics](#)

[SerializeToJson](#)

Serializes a sequence of data model objects into JSON strings.

[SerializeToYaml](#)

Serializes a sequence of data model objects into YAML strings.

[Session](#)

[Sound](#)

[TimeConstrains](#)

Enums

[OptoLEDMode](#)

[OptogeneticsMode](#)

[OptogeneticsRampMode](#)